

<https://doi.org/10.1038/s42004-025-01784-9>

Spitrobot-2 advances time-resolved cryo-trapping crystallography to under 25 ms



Maria Spiliopoulou¹, Caitlin E. Hatton², Martin Kollwe³, Jan-Philipp Leimkohl³, Hendrik Schikora³, Friedjof Tellkamp^{3,4}✉, Pedram Mehrabi^{2,3,4}✉ & Eike C. Schulz^{1,2,3,4}✉

We previously introduced the *spitrobot*, a protein crystal plunging system that enables reaction quenching via cryo-trapping with a time resolution in the millisecond range. Here we present the next generation, *spitrobot-2*, as an integrated benchtop device. User-friendliness has been improved by semi-automatic sample exchange. Moreover, a fully automated shutter shields the liquid nitrogen from the humidified environment, improving sample integrity. Most importantly, the cryo-trapping delay time has been reduced to 23 ms, making *spitrobot-2* twice as fast as the previous generation. This further expands the number of target systems that can be addressed by cryo-trapping time-resolved crystallography. Using 12 crystal structures of three independent model systems, we demonstrate successful cryo-trapping via observation of conformational changes and ligand binding within 25 ms. These improvements increase the convenient access to cryo-trapping, time-resolved X-ray crystallography empowering the MX community with efficient tools to advance research in structural biology.

Time-resolved crystallography (TRX) is a powerful technique to study dynamic events and conformational responses of proteins while they carry out their function. In the past decade, it has undergone a resurgence, primarily triggered by the advent of X-ray free-electron laser (XFEL) sources and the emergence of serial data collection methods^{1–3}. The successful application of serial crystallography at XFELs was quickly adopted at synchrotron facilities, which can provide equivalent data quality for sufficiently large crystals^{4–7}. While serial synchrotron crystallography (SSX) at monochromatic MX beamlines cannot reach the ultra-fast time-scales that are accessible to XFELs, they are, however, the obvious choice when it comes to studying enzymatic mechanisms beyond the microsecond time domain. Their comparably wide distribution goes in hand with a lower competition for beamtime and thus enables a more democratic access, connected with the ability to study a larger number of model systems and develop novel methods that will benefit the field as a whole. The recent developments in time-resolved serial synchrotron crystallography (TR-SSX) have been the subject of a number of reviews^{2,8–14}.

However, in spite of clear fundamental and practical advantages, as well as recent progress in simplifying the method, time-resolved crystallography still requires expert knowledge and thus often remains the niche of experts. Cryo-trapping crystallography is a convenient workaround that permits resolving metastable reaction intermediates by quickly quenching biochemical reactions in liquid nitrogen. Traditionally, slow enzymes with

turnover times in the second's domain or those with repetitive mechanisms have been addressed with cryo-trapping for decades^{15,16}. While true TR-SSX experiments at room- or even physiological temperatures provide a much more complete insight into protein dynamics, and are not affected by vitrification artifacts^{10,17–20}, these can be conveniently complemented by cryo-trapping experiments. Cryo-trapping allows for remote and decoupled sample preparation and data collection, requires fewer crystals for time-resolved studies, which favours challenging systems, e.g., hard to produce or unfavorable crystal size-to-diffraction ratios, and is compatible with low-brightness beamlines or even home sources, thus extending experimental capabilities to a wider range of setups and users. In addition to its conceptual simplicity, a major advantage of cryo-trapping is rooted in its compatibility with established high-throughput infrastructure and automated data-processing routines. However, a common problem in *manual* cryo-trapping lies in its limitation to macroscopic crystals, the comparably large time jitter associated with this, its comparably low reproducibility, especially if time-scales faster than ca. 30 s are aimed for. To address these limitations, we have recently developed the *spitrobot*, an automatic crystal plunger that enables time-resolved crystallography, via cryo-trapping, in the millisecond time domain²¹. Reactions are initiated by the *liquid application method for time-resolved applications* (LAMA), which permits in situ mixing with a minimal amount of substrate solution, while allowing for reaction initiation times in the millisecond time-domain. LAMA nozzles spray picoliter sized droplets

¹University Medical Center Hamburg-Eppendorf (UKE), Hamburg, Germany. ²Institute for Nanostructure and Solid State Physics, Universität Hamburg, Hamburg, Germany. ³Max-Planck-Institute for Structure and Dynamics of Matter, Hamburg, Germany. ⁴These authors jointly supervised this work: Friedjof Tellkamp, Pedram Mehrabi, Eike C. Schulz. ✉e-mail: friedjof.tellkamp@mpsd.mpg.de; pedram.mehrabi@uni-hamburg.de; ec.schulz@uke.de

of ligand solution onto protein crystals, which are exposed to X-rays or cryo-trapped after a pre-defined delay time²². The *spitrobot* is compatible with macroscopic crystals, micro-crystals, and canonical rotation, as well as serial data collection methods. In addition, the *spitrobot* is fully compatible with the SPINE standard (e.g., crystals are directly vitrified inside SPINE pucks), which directly connects cryo-trapping to the high-throughput infrastructure available at most synchrotrons. Another major advantage in cryo-trapping approaches lies in the ability to uncouple sample preparation from data collection, that is, users can prepare their samples well in advance to a beamtime and fully focus on either task.

For the growing user base interested in cryo-trapping crystallography using the *spitrobot*, we have updated the first generation prototype. *Spitrobot-2* is a fully integrated benchtop device with a minimal footprint, conveniently fitting into existing MX-laboratories. Like the first generation, *spitrobot-2* also maintains humidity and temperature conditions during reaction initiation via implementation of an environmental control system, and thereby also permits addressing long delay times²¹. Finally, we have substantially enhanced the user experience by improving the overall design and workflow of *spitrobot-2*. Most importantly, with further hardware improvements, we were able to reduce the minimum delay time by a factor of 2, now enabling cryo-trapping within less than 25 ms.

Results and discussion

An integrated benchtop device

In order to turn the *spitrobot* prototype into a user-friendly device, we aimed to generate a compact setup with minimal clutter that aids the user in the sample preparation process. For an overview incl. abbreviations, please refer to Fig. 1a, b. The footprint of *spitrobot-2* was substantially reduced, now comprising dimensions of $W284 \times H480 \times D316$ mm, with a weight of approximately 15 kg. Two handles on the sides allow for convenient relocation and enable use as a benchtop device. In contrast to the prototype, *spitrobot-2* permits triggering the plunging process from the main device itself, and therefore no longer depends on an external control box. The humidity flow device (HFD), as well as a control computer (for e.g., delay time setup, sample monitoring, automatic image capturing), remains separate.

The front panel comprises the two-hand-control safety switches that were formerly part of the control box (B1, B2). Pressing these buttons triggers reaction initiation and sample plunging according to the settings in the software interface. The two-hand-control solution prevents inadvertent interaction with the sample area during the plunging process and thereby prevents user injury. An additional button (B3) on the right side is used to reset the device in case of an error or as a light switch for sample illumination. To enable high-throughput sample preparation, samples are plunged directly into a SPINE puck. As in the prototype, a modified foam Dewar containing the SPINE puck is placed at the bottom, and its height is adjusted via an integrated LabJack. Sample exchange is executed via a sample exchange dial (SD1) that rotates the SPINE puck inside of the Dewar to access the next sample position (Supplementary Fig. 1). Next to the Dewar, on the lower left of the device, another button (B4) can be used to manually release the pin from the electromagnet on the piston. All buttons (B1–B4) are equipped with a color ring indicating the process or function status, e.g., red for error.

The sample viewing system consists of two cameras, aligned at an angle of 90° to each other. In *spitrobot-2* both cameras can be focused independently via two dials (CD1, CD2) on top of the device. To facilitate non-standard adjustments, such as precise camera centering, the device incorporates a built-in Allen key mechanism, which is automatically engaged and released by a spring-release mechanism, triggered by pressing the button. The LAMA nozzle for reaction initiation is mounted on the left side of *spitrobot-2*, and its precise orientation with respect to the sample position is controlled via three nozzle dials (ND1, ND2, ND3). The nozzle can be locked in a retracted position for sample mounting and quickly released into spraying position by the nozzle release latch. Different nozzle sizes are commercially available, enabling adjustment of the volume of the deposited

substrate solution up to 3 nL/ms (Microdrop LLC, Norderstedt, Germany). The nozzle can be replaced by focusing optics, enabling optical excitation of the sample by a fibre-coupled light source or laser²³. On the right side of *spitrobot-2* the HFD-nozzle is fed through to the sample position. On the back, there are two network ports, four BNC connections (e.g., for the LAMA trigger signal), as well as the main power supply and switch.

A liquid nitrogen level indicator ensures sample quality

Sample quality is key to any TRX experiments, but for cryo-trapping TRX, this also includes the quality of the cryogen. The quality of the vitrified samples critically depends on the steepness of the temperature gradient and the purity of the liquid nitrogen, i.e., its level of water-ice contamination^{24–28}. To support this, integrated temperature sensors warn the user or report an error if the level of liquid nitrogen drops too low, to avoid compromising the vitrification process and affecting the delay time. Liquid nitrogen refilling is possible via a safety funnel, which enables refilling while the Dewar remains in the sample preparation position. To insulate the liquid nitrogen from atmospheric air and the stream of highly humid, warm air from the HFD, *spitrobot-2* is equipped with a liquid nitrogen shutter window that only opens during the plunging period (Fig. 1c), blocking access to the cryogen at all other times. This limits liquid nitrogen evaporation and increases the cryogen life time in the Dewar, but more importantly, it reduces ice contamination, thereby improving sample quality. During the release stage, the window adopts a semi-open position, allowing for swift and reliable sample vitrification.

There is a clear influence of plunge freezing speeds and sample volumes on cooling rates and vitrification outcomes, which provides insights for further optimization. Previous studies on cryofixation emphasize that eliminating the cold gas layer above the liquid nitrogen during plunge freezing enhances cooling rates, thereby improving sample preservation²⁹. Most of the cooling occurs in the cold gas layer above the cryogen during plunge cooling of protein crystals in liquid nitrogen²⁷, and for plunge speeds of 1 ms^{-1} , a 2 cm gas layer is sufficient to dominate the cooling rate^{27,30}. *Spitrobot-2* is capable of efficient plunge freezing at 1.74 ms^{-1} , further improving sample preparation.

Spitrobot-2 enables delay times of under 25 ms

Early vitrification experiments have provided valuable information on the cooling rates required to achieve successful vitrification. For pure water, these need to be as high as 10^6 Ks^{-1} , but this can be reduced by orders of magnitude by the presence of solutes^{24,31}. Initial studies using liquid propane near its melting point as a cryogen achieved a peak cooling rate of about $370,000 \text{ Ks}^{-1}$, determined by 25–75 μm thermocouples immersed at a rate of 2 ms^{-1} ³². Further research has highlighted significant differences in the cooling rates of different cryogenic media. In particular, the cooling rates in boiling liquid nitrogen have been reported to be about 50 times lower than those in, for example, 90 K ethane³³, allowing vitrification without cryo-protectant. Although liquid nitrogen requires the use of cryo-protectants, the reduced safety concerns compared to flammable cryogenics have made it the primary cryogen for X-ray crystallography^{25,34} and thus the primary choice for *spitrobot-2*.

The electro-pneumatic piston of the prototype has been replaced by a motorised solution, which in *spitrobot-2* drives the sample into a vial submerged in liquid nitrogen. In total, it is equipped with two linear motors (Faulhaber LM series), the first of which is used for the plunging process, while the second controls a shutter between the dispensing chamber and the Dewar. The plunger drives (Fig. 1c) the samples into the cryogen at a speed of 1.74 ms^{-1} (40 mm in 23 ms). To experimentally validate the minimum delay time that can be achieved with *spitrobot-2*, we measured the temperature evolution by plunging a 13 μm thermocouple into liquid nitrogen, with and without applying the humidity stream, as described previously²¹.

The metallic thermocouple sensor (nickel, type K) used for the cooling rate measurement has a diameter of 13 μm , while the crystals used in this study have diameters of around 20 μm . Based on their respective specific heat capacities (Table 1) we have determined that the equilibration time for

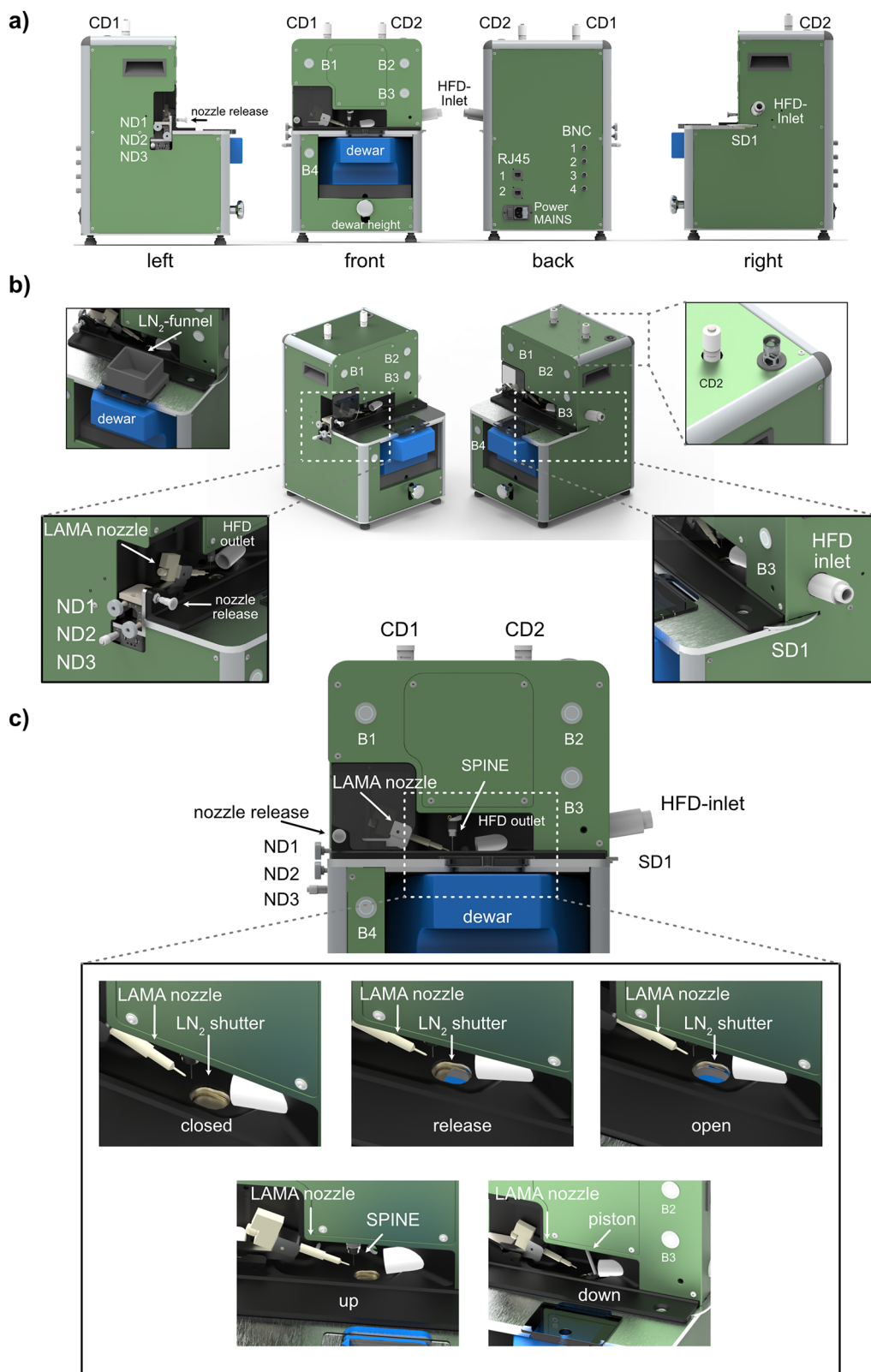


Fig. 1 | The redesigned *spitrobot-2*. **a** All side views of *spitrobot-2*. **b** 2-sided isometric view of *spitrobot-2* with closeups of the liquid nitrogen (LN₂) funnel, the allen key lift, the nozzle dials (ND1-ND3), and the sample dial (SD). **c** The front-side of

spitrobot-2, with closeups of the different liquid nitrogen shutter steps and the two main positions of the piston. Abbreviations: CD camera dial, ND nozzle dial, B button, SD sample dial, HFD humidity flow device.

the thermocouple as well as for protein crystals is in the micro-second domain, well below the experimentally measured quench time (Materials and Methods). Therefore, the thermocouple derived cooling times are an appropriate estimate for the protein crystals used in this study.

As these data clearly show, *spitrobot-2* allows complete vitrification of the samples within 23 ms, which is about twice as fast as the prototype (Fig. 2a). For crystal sizes that are typically still useful for diffraction experiments at synchrotron sources (about 10–20 μm), this delay time now

Table 1 | Thermodynamic parameters comparing the thermocouple and protein crystals

	thermocouple (nickel)	protein crystal (lysozyme)	
diameter of the sphere d (μm)	13	20	40
specific heat capacity c_p ($\frac{\text{J}}{\text{kg}\cdot\text{K}}$)	444 ⁷³	1750 ^{58,74,75}	
density ρ ($\frac{\text{kg}}{\text{m}^3}$)	8908 ⁷⁶	1240 ^{77,78}	
thermal conductivity κ (W/mK)	90 ⁷⁹	0.42 ⁸⁰	
heat capacity of the sphere C ($\frac{\text{J}}{\text{K}}$)	4.5	9.1	72.7
thermal diffusivity α (m^2/s)	$23 \cdot 10^{-6}$	$0.194 \cdot 10^{-6}$	
equilibration time t_e (μs)	0.55	155	620

approaches the theoretical lower diffusion time limit for fast diffusing ligands such as glucose (15 ms)^{10,22,35}. For practical purposes, it is important to consider that the ligand diffusion rate into protein crystals does not only depend on crystal packing and solvent channels, but also on the viscosity of the cryoprotectant. All these factors contribute to the adjustment of solvent viscosity and hydration levels that affect diffusion kinetics and ligand accessibility within the crystal lattice³⁶. Notably, by using the HFD in *spitrobot-2* the cooling time measured by the thermocouple could be further reduced to 3.5 ms (Fig. 2 b). This equates to a cooling rate of approximately 57,000 Ks⁻¹. While complete cooling to 77 K occurs on the order of a few milliseconds, the functionally relevant 'quench' point is when the crystal passes below approximately 200 K, at which point enzymatic activity ceases^{17,37–39}. For 20 μm crystals, this occurs within ca. 1–2 ms, and ca. 1 ms to reach the glass transition temperature, consistent with a recent study⁴⁰. We anticipate that by selecting alternative cryogenic media, higher cooling rates may be possible with *spitrobot-2* in the future, potentially bypassing the need for cryo-protection.

A recent article describes instrumentation similar to *spitrobot-2*⁴⁰. However, a side-by-side comparison of our system with the *mix-and-quench* approach of Indergaard et al. reveals some substantive differences:

Spitrobot-2 achieves a minimum delay time of 23 ms, comprising the time from ligand application to complete vitrification. On the other hand, the *mix-and-quench* approach claims 8 ms delays by referring to the interval between ligand addition at the loop and immersion in liquid nitrogen, but excluding the vitrification time. Furthermore, although an 8 ms delay time has been reported, this has unfortunately not been validated by structural data, as there is no evidence of ligand binding below 150 ms⁴⁰.

Moving towards the same objective and acknowledging the current needs in the field for more efficient time-resolved crystallography tools, *spitrobot-2* was developed to overcome some of the limitations observed in the previous generation (*spitrobot-1*), particularly regarding time delay. As described previously²¹, *spitrobot-1* is equipped with an air valve, and the motion of the piston is therefore delayed by 20 ms. *Spitrobot-2* is instead equipped with an electrical linear motor, enabling faster plunging of the sample to the liquid nitrogen (Fig. 3). This enables cryo-trapping experiments at delay times close to the theoretical diffusion limits.

Sample preparation workflow

The original workflow of the *spitrobot* has been redesigned to provide a more user-friendly, streamlined experience that speeds up the entire sample preparation process. To start a sample preparation, the pre-cooled Dewar containing the SPINE puck is now placed on an integrated lifting mechanism. Once the Dewar has been inserted, precise positioning within the *spitrobot-2* is aided by a bright LED guiding beam which helps the user to align the puck with respect to the axis of the instrument's rotational mechanism (Supplementary Fig. 1). Once in position, the Dewar is lifted by the labjack until it seals against the insulated and actively heated plate separating the Dewar from the upper part of the machine. This has two

advantages: firstly, it helps to protect the operator during sample vitrification; secondly, it also greatly reduces ice formation within the liquid nitrogen as the Dewar is now isolated from both the atmosphere and the moisture stream. Thus, with these improvements the use time is approximately 30 min before refilling of liquid nitrogen is required. However, it is good common practice to top-up the liquid nitrogen after 3–5 samples have been plunged, to reduce sample variability. Sample access to the liquid nitrogen is only possible via the integrated shutter (*see above*). Next, the semi-manual sample dial on the right-hand side of the *spitrobot-2* is used to position the start vial in the immersion position. The alignment of the dispenser nozzle is supported by an integrated dual camera system (IDS imaging). Now the desired delay time can be set in the control software and the system is ready to start the plunging process by pressing the two-hand-control safety switches. Whilst retaining the two-handed safety switch solution implemented in the prototype, these have been integrated into the main unit so that *spitrobot-2* is no longer dependent on an external control box. Moreover, colour-coded signals integrated into the control buttons provide the user with immediate feedback about the plunging process, e.g., green and blue to let the operator know when to press or release buttons, or a white flash at the moment of sample excitation.

Ligand binding events and conformational changes

To demonstrate the versatile applicability of *spitrobot-2* we demonstrate reaction initiation and fast delay times in three independent model systems, obtained by single crystal rotation data collection: *Streptomyces rubiginosus* xylose isomerase (XI), human insulin (HI), and bacteriophage T4 lysozyme (T4L).

Xylose isomerase ligand binding trapped in 25 ms. To allow direct comparison to the first generation of *spitrobot*, we used xylose isomerase (XI) as a model system. XI is of great commercial interest as it catalyses the reversible interconversion of α -D-glucose to β -D-fructose^{41,42}. We have previously shown that glucose can bind to XI within 15 ms, which approximates theoretical diffusion times for crystals of this size^{10,22,35}. To initiate substrate binding, we sprayed a 2 M glucose solution onto 20 μm XI crystals loaded on a micromesh before plunging them into liquid nitrogen with delay times of 25 ms and 50 ms, respectively. The difference density map in the active site shows that with *spitrobot-2* we can capture the diffusion of glucose after the minimum delay time of 25 ms (Fig. 4). This means that our current data confirm previous observations on glucose diffusion kinetics in protein crystals²². The onset of glucose binding was reported to occur at 15 ms using the same enzymatic system via SSX, and is consistent with the identification of a partially occupied state at 25 ms. In both studies, the ligand is observed in the active site, in the same conformation.

A pH jump triggers conformational changes in human insulin.

Human insulin (HI) is a 5.8 kDa peptide hormone produced by β -pancreatic cells, which facilitates the absorption of glucose from the blood into tissues. It was one of the first proteins to be isolated and crystallographically studied^{143–46}. In its active form, HI consists of 51 amino acids divided into two polypeptide chains: A and B (21 and 30 amino acids, respectively). The secondary structure of insulin shows two almost antiparallel α -helices in chain A and a α -helix followed by a turn and a β -strand in chain B⁴⁷. After a century of crystallographic research, insulin has been identified as one of the most polymorphic proteins, displaying a range of molecular conformations, crystal forms, and binding capacities for metals and ligands⁴⁸. Several cubic insulin crystal structures have been analysed under different conditions, such as crystals in 0.1 M sodium salt solutions with pH values ranging from 7.0 to 10.0⁴⁹ and in 1 M Na₂SO₄ at pH 5.0 to 11.0⁵⁰. It has previously been reported that cubic insulin crystals undergo a structural transition when the pH is lowered from 9.0 to 5.0⁵⁰, while their crystal packing is not affected. The conformational change of glutamate at position 13 of chain B (GluB13) is accompanied by the binding of sulphate ions close to phenylalanine at position 1 of chain B

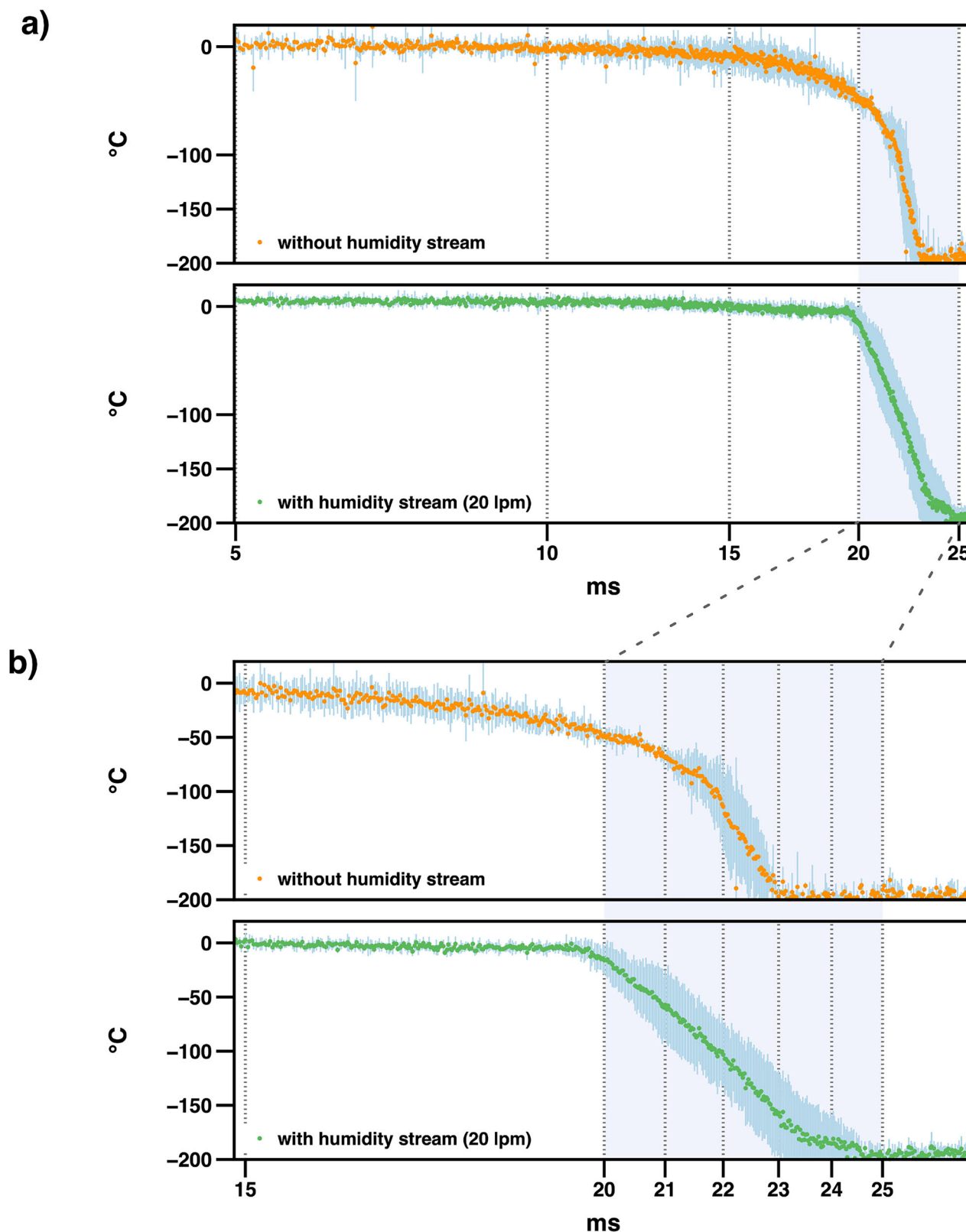


Fig. 2 | Average *spitrobot-2* vitrification time in liquid nitrogen. The temperature profile is shown without (orange, top-panel) and with (green, bottom panel) the humidity flow device, respectively. The standard deviation is shown as blue lines. Overall, the plunge process took 22.6 ± 0.4 and 23.6 ± 0.7 ms, respectively. **a** A global view of the average temperature evolution between 5 and 25 ms. Irrespective of the humidity flow device, the temperature drops to below -196°C , within less than

25 ms. **b** A closeup of the vitrification time between 15 and 25 ms demonstrates that by application of the humidity flow device, the temperature remains stable until the sample enters the liquid nitrogen. The vitrification process in liquid nitrogen itself, then takes around 3.5 ms and thus is about $2 \times$ faster than the vitrification time previously reported for the first generation.

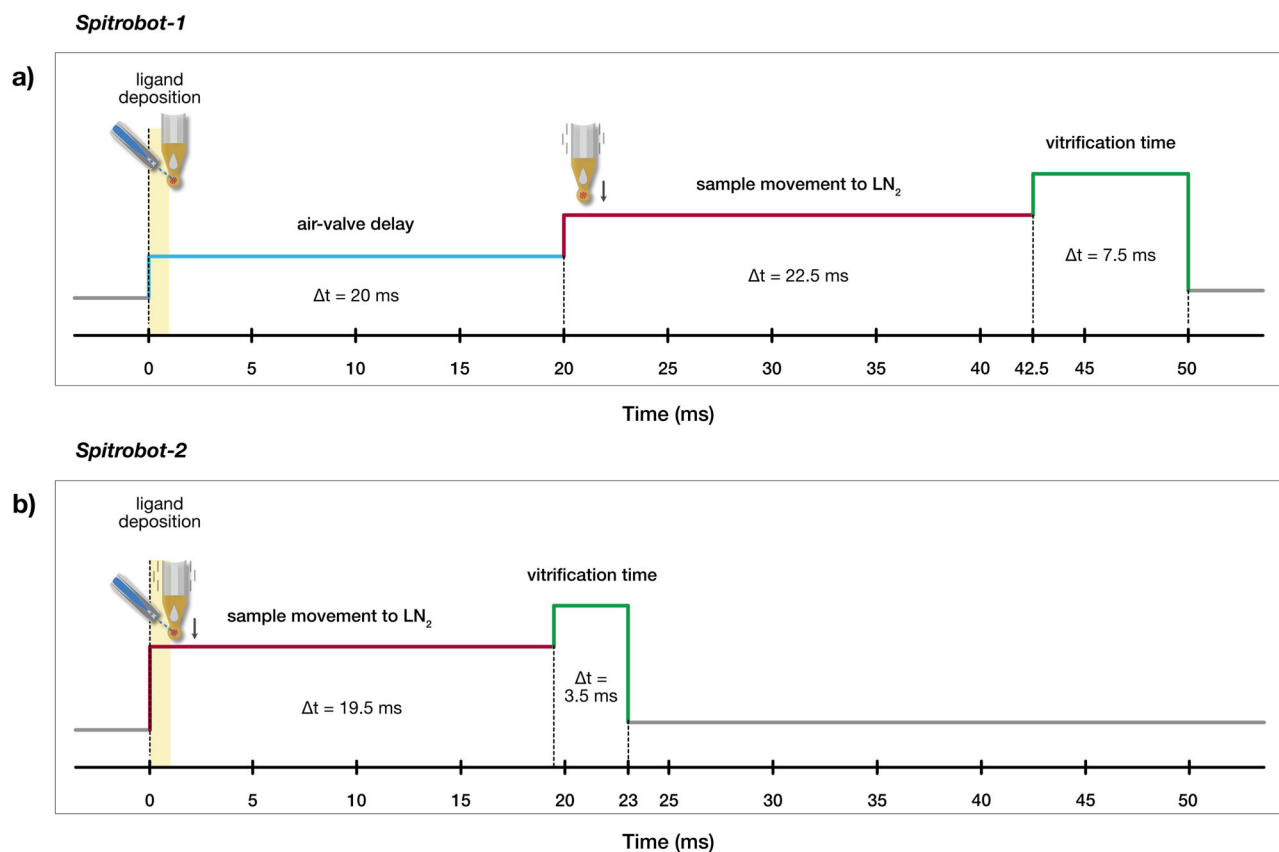


Fig. 3 | Comparison of the delay time between the spitrobot generations. **a** *Spitrobot-1*: The minimum total delay time is 50 ms: air-valve delay (20 ms), mechanical piston motion (22.5 ms), and vitrification time (7.5 ms); the ligand

deposition (yellow) occurs during the air-valve delay. **b** *Spitrobot-2*: The minimum total delay time is 23 ms: electrical motor (plunger) motion (19.5 ms) and vitrification time (3.5 ms); the ligand deposition (yellow) occurs during plunger motion.

(PheB1). Here, we aimed to induce this effect by a pH jump with *spitrobot-2*.

To this end, we sprayed a low pH sulphate buffer (1 M CH_3COONa , pH 4.5, 1 M Na_2SO_4 , 15% (v/v) ethylene glycol) onto the insulin crystals and followed the motion of GluB13 and SO_4 binding as a function of time. As can be clearly seen from the difference density maps, both the structural transition and the SO_4 binding can be observed at delay times as short as 25 ms (Fig. 5; for refined occupancies for GluB13 and SO_4 see Supplementary Table 2).

Backbone changes in bacteriophage T4 lysozyme (T4L). The biological function of T4L is to hydrolyse the β -1,4 linkage between N-acetylmuramic acid and N-acetylglucosamine in bacterial peptidoglycan to aid cell lysis^{51,52}. The L99A mutation of T4L forms an internal non-polar cavity of about 150 Å³, which allows the binding of various ligands after large conformational changes that open T4L to allow ligand binding^{53,54}. As a result, this particular mutant serves as a model system for studying ligand binding at allosteric sites. Previous studies focusing on benzene binding pathways have proposed that structural rearrangements of the cavity surrounding the helices allow the ligand to enter the binding site^{53–55}. To explore whether we could induce and capture larger conformational changes, we used indole as a ligand and observed structural changes after 1 s and 10 s, respectively. Our *spitrobot-2* cryo-trapping data clearly show that structural snapshots of the binding event can be obtained at delay times (1 s) inaccessible by manual methods (Fig. 6). In addition to indole binding, we clearly observe a conformational change in the backbone of the F-helix of T4L in the indole-bound states, which is shifted by about 1.8 Å compared to the apo state (Fig. 6e).

Conclusions

Spitrobot-2 is an integrated benchtop solution for cryo-trapping protein crystals with a delay time of less than 25 ms. It is suitable for monitoring ligand binding events as well as conformational changes of main and side chains. As the average enzyme has a turnover number of around 10 s^{−156}, a large number of systems should be amenable to cryo-trapping TRX, e.g., to track catalytic intermediates that are inaccessible by conventional methods. This is underlined by the fact that six highly different enzymes have now been studied by both generations of the *spitrobot*^{21,23}—testament to the device's versatility and reliability. Combined with its small footprint and conceptual simplicity, these features should make *spitrobot-2* an attractive asset for MX labs wishing to standardize their vitrification process or explore cryo-trapping TRX. Its strength is further emphasized by its full compatibility with the high-throughput workflows found on most MX beamlines.

Methods

Temperature measurements

To assess the temperature evolution during the plunging process, we attached a 13 μm type K thermocouple (KFT-13-200-200, As One Corporation) to a canonical SPINE pin. The thermocouple was mounted on a SPINE pin at a distance of 22.0 mm from the base, within ±1 mm at the sample position. This modified test sample was plunged into liquid nitrogen 3 times without, and 5 times with the humidity flow device, respectively. Temperature corresponding voltage measurements were recorded at 50 μs intervals. Furthermore, in order to reduce noise, we used a 100 μs low pass filter to connect the thermocouple to the oscilloscope (bandwidth: 350 MHz). Consequently, the readout lag is predominantly influenced by the low pass filter, but remains negligible in comparison to the overall

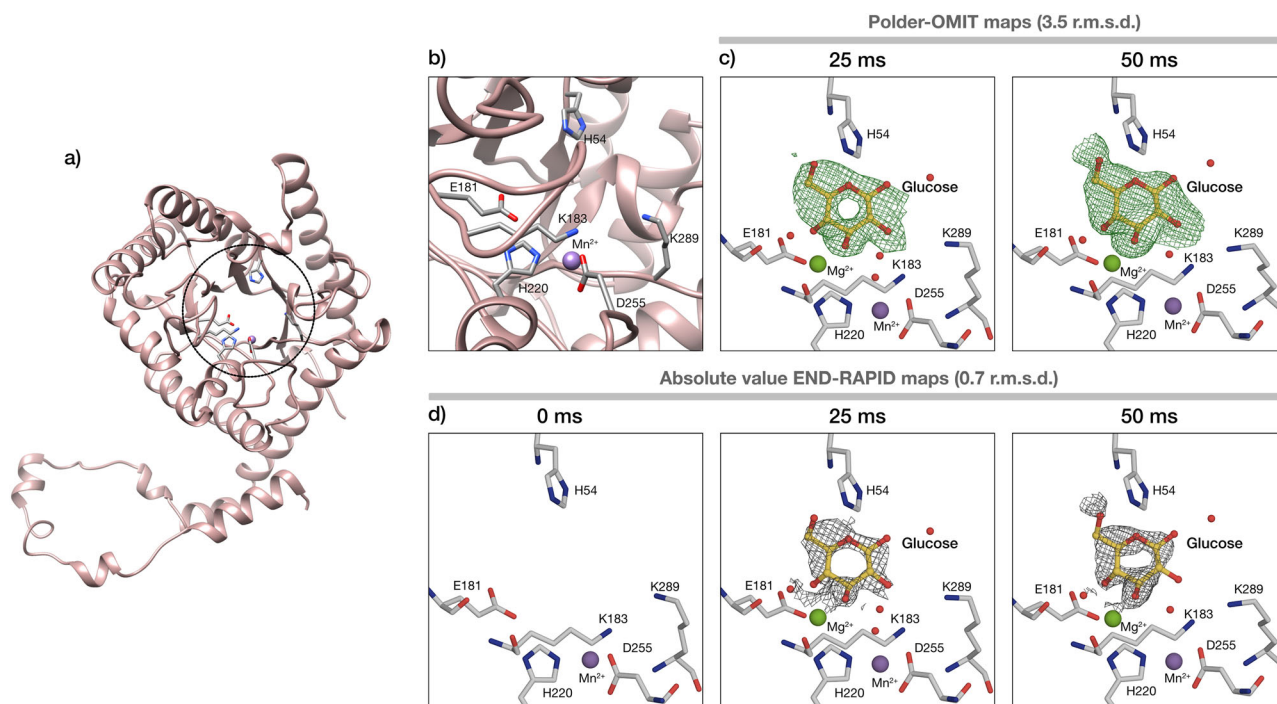


Fig. 4 | Glucose binding in the active site of XI. **a** Cartoon representation of the XI structure; the active site is indicated by a circle. **b** Closeup of the active site and surrounding residues. **c, d** Polder-OMIT maps and absolute value END-RAPID

maps for glucose binding after 25 and 50 ms of soaking, respectively. Mg²⁺ and Mn²⁺ are represented as green and purple spheres correspondingly.

vitrification time. In this context, vitrification time refers to the experimentally measure cooling time of the thermocouple.

Equilibration times. To clarify whether the thermocouple provides a sufficiently good approximation for the cooling of protein crystals, we determined the equilibration time of the thermocouple and the equilibration time of lysozyme crystals of different sizes. We consider a homogeneous spherical sample with thermal diffusivity α . The equilibration time t_e is defined as the time required for the temperature at the center of the sphere to reach 90% of the imposed temperature change at its boundary. Therefore,

$$t_e \approx 0.3 \cdot \frac{r^2}{\alpha},$$

with r the radius of the sphere⁵⁷. The thermal diffusivity α is given by the ratio of the thermal conductivity κ to the density of the material ρ and its specific heat capacity c_p .

$$\alpha = \frac{\kappa}{\rho \cdot c_p},$$

The heat capacity C for the samples is given by

$$C = V \cdot \rho \cdot c_p,$$

wherein V is the volume of a sphere ($\frac{4}{3}\pi r^3$), r is the radius of the sphere, ρ is the density of the material, and c_p is the specific heat capacity. Based on the parameters in Table 1, the equilibration time within the 13 μm thermocouple is 0.55 μs , compared to 155 μs for 20 μm lysozyme crystals.

It is important to note that this estimate reflects the internal equilibration time within the sample, not the total quench time. This is set by the heat transfer from the surface of the sphere into the cryogen. Heat transfer from protein crystals during plunge-cooling into liquid nitrogen has been addressed in detail by Kriminski et al. They worked out that for sufficiently small samples (<500 μm), heat transfer is limited by the liquid

boundary layer⁵⁸. This is in agreement with earlier work, showing that internal heat conduction timescales are significantly shorter, thereby justifying the dominance of surface heat transfer in controlling the cooling dynamics⁵⁷.

While the thermal diffusivity of lysozyme crystals is comparatively low with respect to the nickel thermocouple, the internal heat conduction is still well below the experimentally measured quench time. Therefore, the experimentally determined cooling of the nickel thermocouple provides a reliable benchmark for the effective quench rate using our plunge conditions.

Protein purification and crystallization

Xylose isomerase. *Streptomyces rubiginosus* xylose isomerase in pET24a vector (Genscript), was expressed in BL21 DE3 E. coli cells, purified by affinity (Ni) and size exclusion chromatography with a final buffer of 0.05 M Tris pH 8.5, 0.15 M NaCl, with crystallization as reported in previous studies²². Briefly: 80 mg mL⁻¹ XI was combined with an equal volume of crystallization buffer (35% (w/v) PEG 3350, 0.2 M LiSO₄, and 0.01 M HEPES/NaOH, pH 7.5) via SpeedVac crystallization as previously described⁵⁹.

Human insulin. Human insulin purchased from Roche (LOT: 70272900) was diluted in 0.05 M Na₂HPO₄, 0.001 M Na₂EDTA, pH 11.0 (adjusted) to a final concentration of 30 mg mL⁻¹. Batch crystallization was performed in a total volume of 160 μL by mixing equal volumes of HI and crystallization buffer consisting of 1.0 M NaH₂PO₄/K₂HPO₄, pH 9.0⁶⁰.

Bacteriophage T4 lysozyme. T4-L99A mutant lysozyme gene was cloned in pET29b(+) vector (Genscript) and expressed in BL21 DE3 E. coli cells. The protein was purified by ion exchange and size exclusion chromatography with a final buffer of 0.05 M NaH₂PO₄/Na₂HPO₄, pH 5.5, 0.1 M NaCl, 0.002 M EDTA⁶¹. Batch crystallization was performed in a total volume of 20 μL by mixing equal volumes of protein at 22 mg mL⁻¹ and crystallization buffer consisting of 4.0 M NaH₂PO₄/K₂HPO₄, pH 7.0, 0.1 M 1,6-hexanediol, and 0.15 M NaCl (PDB ID: 3K2R).

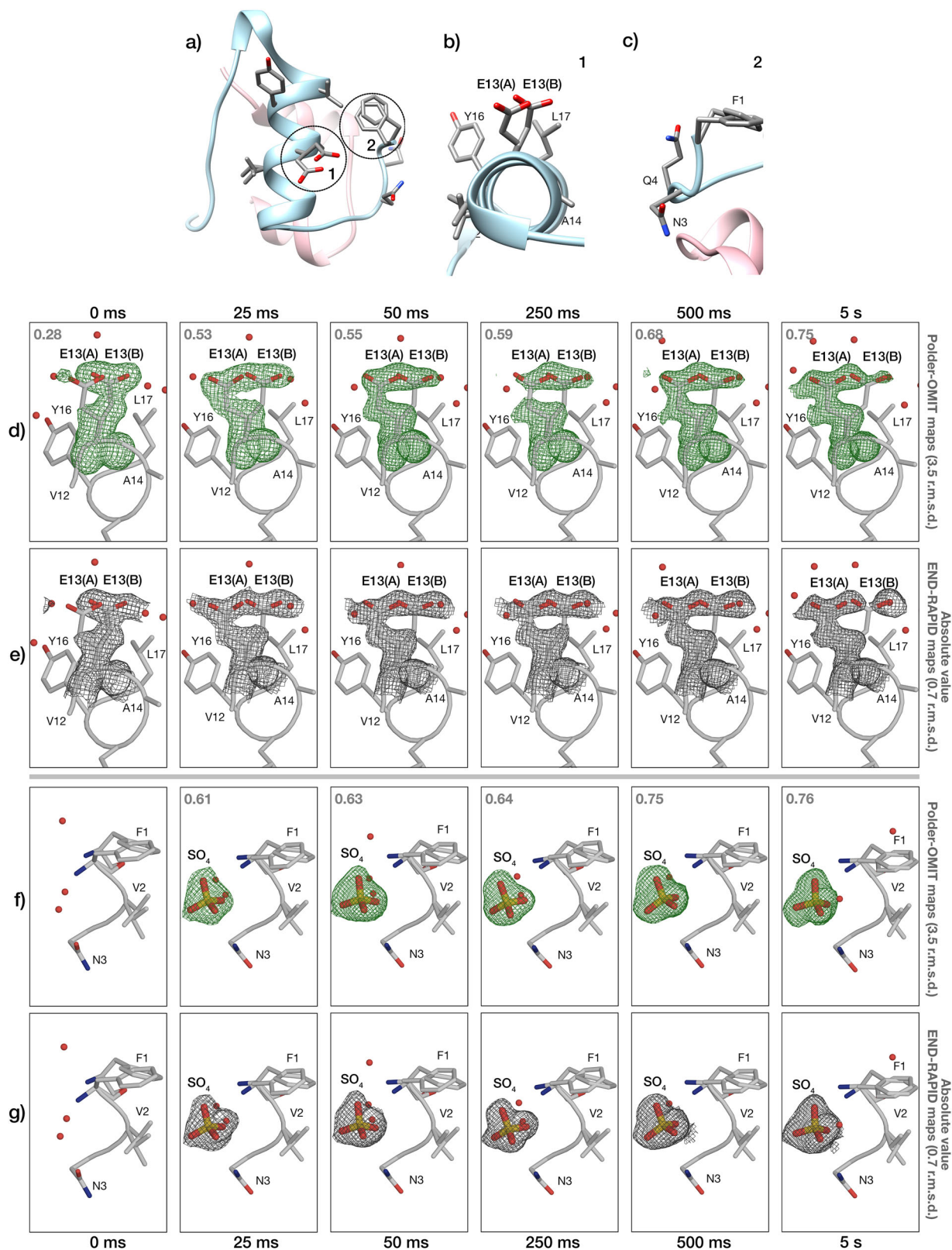


Fig. 5 | pH-jump effect on HI. **a** Cartoon representation of HI structure; the two chains are illustrated in pink (chain A) and blue (chain B), while the residues in the GluB13 region (1) and the SO_4 binding region (2) are represented in sticks. **b** Closeup of GluB13 and surrounding residues for HI at pH 9.0. **c** Closeup of the SO_4 binding region at pH 9.0 (apo state). **d, e** Polder-OMIT maps and absolute value END-

RAPID maps for the effect of pH 4.5 on the GluB13 after 25 ms to 5 s; the occupancies of conformation A are indicated in grey. **f, g** Polder-OMIT maps and absolute value END-RAPID maps for the effect of pH 4.5 on SO_4 binding after 25 ms to 5 s soaking; the occupancies for SO_4 are indicated in grey.

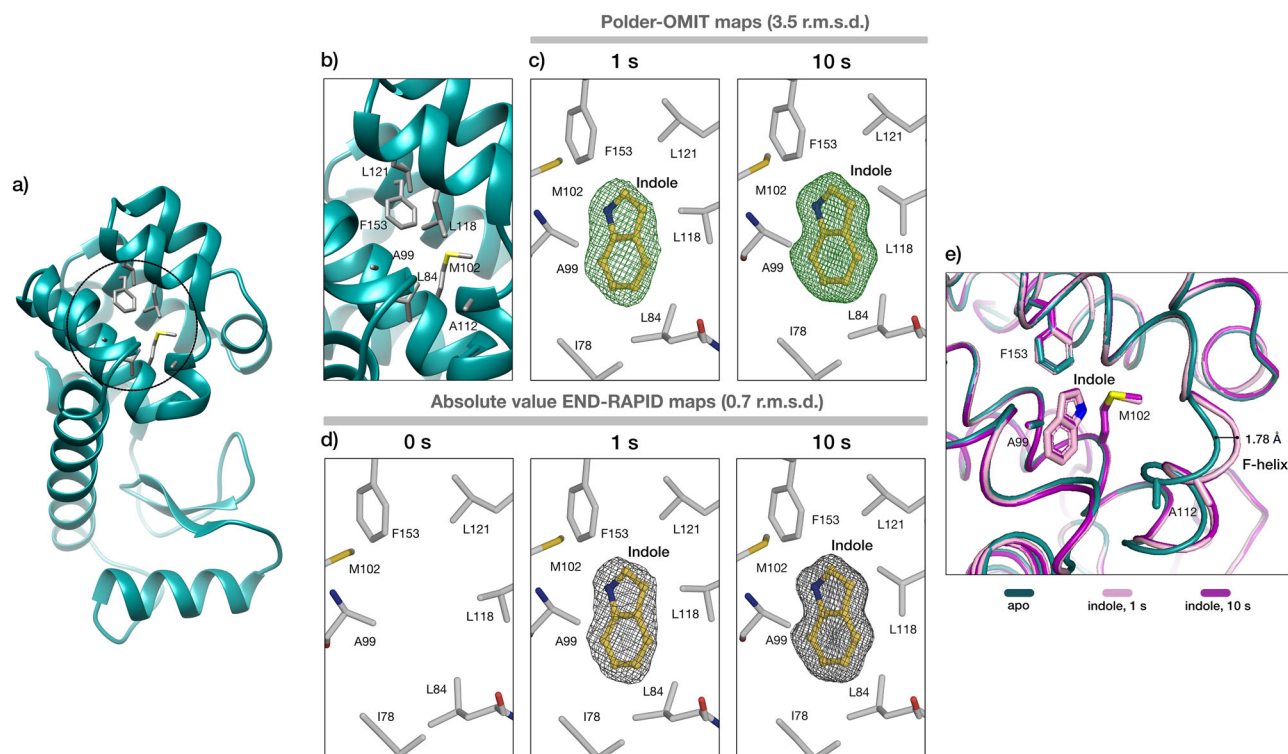


Fig. 6 | Indole binding in the non-polar cavity of T4L-L99A. **a** Cartoon representation of T4L-L99A structure; the cavity is indicated by a circle. **b** Closeup of the cavity and the surrounding residues. **c, d** Polder-OMIT maps and absolute value END-RAPID maps for indole binding after 1 s and 10 s of soaking time. **e** T4L-L99A

cavity closeup. Ribbon representation of protein cavity in the absence (teal) and presence of indole after 1 s (pink) and 10 s (magenta) soaking time. The conformational change of the *F*-helix can be observed for the bound states of T4L, leading to a further opening of the cavity (1.78 Å) in the presence of the ligand.

Table 2 | *Spitrobot-2* samples and environmental details

Sample	XI	HI	T4L-L99A
Crystal size (μm)	20	20	20
Temperature (°C)	25	25	25
Relative humidity (%)	95	95	95
Ligand	Glucose	Na ₂ SO ₄	Indole
LAMA solution	0.01 M HEPES pH 7.5, 0.1 M LiSO ₄ , 0.1 M MgCl ₂ , 2 M glucose	1 M CH ₃ COONa pH 4.5, 1 M Na ₂ SO ₄ , 15% (v/v) ethylene glycol	2.0 M NaH ₂ PO ₄ /K ₂ HPO ₄ , pH 7.0, 0.15 M NaCl, 0.1 M 1,6-hexanediol, 6% (v/v) DMSO, 0.05 M indole, 20% (v/v) glycerol

***Spitrobot-2* experiment parameters**

Crystals of the three individual protein systems were used for the *spitrobot-2* experiments. The conditions and solutions used in each case are mentioned in Table 2.

X-ray data collection

Single crystal rotation datasets were collected at beamline P14 of PETRA-III (DESY, Hamburg) at the EMBL Hamburg outstation as well as at beamline ID30A-3 of ESRF (Grenoble, France). Data were collected on an Eiger2 CdTe 16 M detector at 12.7 keV and an Eiger1 X 4 M detector at 12.81 keV, respectively for the aforementioned beamlines. Data collection parameters for each individual protein case are summarized in Supplementary Tables 1, 2 and 3 for XI, HI, and T4L-L99A, correspondingly.

Data processing

Diffraction data were automatically processed with autoPROC using StarAniso^{62–65}. Molecular replacement was conducted in Phaser using PDB-ID 8AWS, 1B18, and 4W51, respectively, as starting models for XI, HI, and T4L-L99A⁶⁶. Structures were refined using iterative cycles of REFMAC or phenix.refine and coot^{67–69}. Molecular images were generated in PyMol⁷⁰.

Different map types were generated in order to examine lower occupancy features; Polder maps were generated to assess weak or ambiguous electron density by omitting both the region of interest and surrounding bulk solvent during map calculation using phenix.polder⁷¹. This approach enhances the visibility of low-occupancy features, such as ligands and alternative conformations, by reducing bulk-solvent model bias and preventing artificial density flattening. END-RAPID maps estimate local noise directly from experimental and model-derived errors⁷². These maps are presented on an absolute scale and visualized using a redefined noise threshold, enabling detection of weak but genuine electron density that may be obscured in conventional maps, and enables direct comparison of electron density features, since maps are on an absolute scale.

Reporting summary

Further information on research design is available in the Nature Portfolio Reporting Summary linked to this article.

Data availability

The data that support this study are available from the corresponding authors upon request. All crystallographic data have been deposited in the

Protein Data Bank (PDB) under the following accession numbers: 9R45, 9R46, 9R47, 9R48, 9R49, 9R4A, 9R4B, 9R4C, 9R4E, 9R4F, 9R4G, 9R4H. Further details are available in Supplementary Tables 1–3. Data related to temperature measurements for the determination of spitrobot-2 vitrification time are also available (Supplementary Data 1).

Due to a patent, detailed designs or a bill of materials cannot be made available at this point. The device is currently available on a limited basis for academic use. Researchers interested in obtaining a unit or collaborating on its use are invited to contact the corresponding authors. We are also evaluating options for broader distribution, including through a potential startup or dedicated production for interested groups.

Received: 3 June 2025; Accepted: 30 October 2025;

Published online: 20 November 2025

References

- Chapman, H. N. X-ray free-electron lasers for the structure and dynamics of macromolecules. *Annu. Rev. Biochem.* **88**, 35–58 (2019).
- Wilson, M. A. Mapping enzyme landscapes by time-resolved crystallography with synchrotron and x-ray free electron laser light. *Annu. Rev. Biophys.* **51**, 79–98 (2022).
- Barends, T. R., Stauch, B., Cherezov, V. & Schlichting, I. Serial femtosecond crystallography. *Nat. Rev. Methods Prim.* **2**, 59 (2022).
- Gati, C. et al. Serial crystallography on in vivo grown microcrystals using synchrotron radiation. *IUCrJ* **1**, 87–94 (2014).
- Stellato, F. et al. Room-temperature macromolecular serial crystallography using synchrotron radiation. *IUCrJ* **1**, 204–212 (2014).
- Weinert, T. et al. Serial millisecond crystallography for routine room-temperature structure determination at synchrotrons. *Nat. Commun.* **8**, 542 (2017).
- Mehrabi, P. et al. Serial femtosecond and serial synchrotron crystallography can yield data of equivalent quality: a systematic comparison. *Sci. Adv.* **7**, eabf1380 (2021).
- Pearson, A. R. & Mehrabi, P. Serial synchrotron crystallography for time-resolved structural biology. *Curr. Opin. Struct. Biol.* **65**, 168–174 (2020).
- Hough, M. A. & Owen, R. L. Serial synchrotron and XFEL crystallography for studies of metalloprotein catalysis. *Curr. Opin. Struct. Biol.* **71**, 232–238 (2021).
- Schulz, E. C., Yorke, B. A., Pearson, A. R. & Mehrabi, P. Best practices for time-resolved serial synchrotron crystallography. *Acta Crystallogr. Sect. D Struct. Biol.* **78**, 14–29 (2022).
- Malla, T. N. & Schmidt, M. Transient state measurements on proteins by time-resolved crystallography. *Curr. Opin. Struct. Biol.* **74**, 102376 (2022).
- Nam, K. H. Guide to serial synchrotron crystallography. *Curr. Res. Struct. Biol.* **7**, 100131 (2024).
- Henkel, A. & Oberthür, D. A snapshot love story: What serial crystallography has done and will do for us. *Biol. Crystallogr.* **80**, 563–579 (2024).
- Banari, A. et al. Advancing time-resolved structural biology: latest strategies in cryo-EM and X-ray crystallography. *Nat. Methods* **22**, 1–16 (2025).
- Hekstra, D. R. Emerging time-resolved x-ray diffraction approaches for protein dynamics. *Annu. Rev. Biophys.* **52**, 255–274 (2023).
- Moffat, K. & Henderson, R. Freeze trapping of reaction intermediates. *Curr. Opin. Struct. Biol.* **5**, 656–663 (1995).
- Halle, B. Biomolecular cryocrystallography: structural changes during flash-cooling. *Proc. Natl. Acad. Sci. USA* **101**, 4793–4798 (2004).
- Fraser, J. S. et al. Accessing protein conformational ensembles using room-temperature X-ray crystallography. *Proc. Natl. Acad. Sci. USA* **108**, 16247–16252 (2011).
- Bock, L. V. & Grubmüller, H. Effects of cryo-EM cooling on structural ensembles. *Biophys. J.* **121**, 148a (2022).
- Schulz, E. C. et al. Probing the modulation of enzyme kinetics by multi-temperature, time-resolved serial crystallography. *Nat. Commun.* **16**, 6553 (2025).
- Mehrabi, P. et al. Millisecond cryo-trapping by the spitrobot crystal plunger simplifies time-resolved crystallography. *Nat. Commun.* **14**, 2365 (2023).
- Mehrabi, P. et al. Liquid application method for time-resolved analyses by serial synchrotron crystallography. *Nat. methods* **16**, 979–982 (2019).
- Kondo, Y. et al. Apo-state structure of the metabotropic glutamate receptor 5 transmembrane domain obtained using a photoswitchable ligand. *Protein Sci.* **34**, e70104 (2025).
- Dubochet, J. et al. Cryo-electron microscopy of vitrified specimens. *Q. Rev. Biophys.* **21**, 129–228 (1988).
- Garman, E. F. & Schneider, T. R. Macromolecular cryocrystallography. *Appl. Crystallogr.* **30**, 211–237 (1997).
- Kriminski, S., Caylor, C., Nonato, M., Finkelstein, K. & Thorne, R. Flash-cooling and annealing of protein crystals. *Biol. Crystallogr.* **58**, 459–471 (2002).
- Warkentin, M., Berejnov, V., Hussein, N. S. & Thorne, R. E. Hyperquenching for protein cryocrystallography. *J. Appl. Crystallogr.* **39**, 805–811 (2006).
- Clinger, J. A., Moreau, D. W., McLeod, M. J., Holyoak, T. & Thorne, R. E. Millisecond mix-and-quench crystallography (MMQX) enables time-resolved studies of PEPCK with remote data collection. *IUCrJ* **8**, 784–792 (2021).
- Ryan, K. P. Cryofixation of tissues for electron microscopy: a review of plunge cooling methods. *Scanning Microsc.* **6**, 8 (1992).
- Berejnov, V., Hussein, N. S., Alsaied, O. A. & Thorne, R. E. Effects of cryoprotectant concentration and cooling rate on vitrification of aqueous solutions. *J. Appl. Crystallogr.* **39**, 244–251 (2006).
- Warkentin, M., Sethna, J. P. & Thorne, R. E. Critical droplet theory explains the glass formability of aqueous solutions. *Phys. Rev. Lett.* **110**, 015703 (2013).
- Costello, M. J. Cryo-electron microscopy of biological samples. *Ultrastruct. Pathol.* **30**, 361–371 (2006).
- Ravelli, R. B. G. et al. Cryo-EM structures from sub-nl volumes using pin-printing and jet vitrification. *Nat. Commun.* **11**, 2563 (2020).
- Pflugrath, J. Practical macromolecular cryocrystallography. *Struct. Biol. Crystall. Commun.* **71**, 622–642 (2015).
- Schmidt, M. et al. Mix and inject: reaction initiation by diffusion for time-resolved macromolecular crystallography. *Adv. Condensed Matter Phys.* **2013**, 167276 (2013).
- Geremia, S., Campagnolo, M., Demitri, N. & Johnson, L. N. Simulation of diffusion time of small molecules in protein crystals. *Structure* **14**, 393–400 (2006).
- Ringe, D. & Petsko, G. A. The ‘glass transition’ in protein dynamics: What it is, why it occurs, and how to exploit it. *Biophys. Chem.* **105**, 667–680 (2003).
- Frauenfelder, H. et al. A unified model of protein dynamics. *Proc. Natl. Acad. Sci. USA* **106**, 5129–5134 (2009).
- Doster, W. The protein-solvent glass transition. *Biochim. Biophys. Acta (BBA) Proteins Proteom.* **1804**, 3–14 (2010).
- Indergaard, J. A. et al. Instrumentation and methods for efficient time-resolved X-ray crystallography of biomolecular systems with sub-10 ms time resolution. *IUCrJ* **12**, 372–383 (2025).
- Marshall, R. O. & Kooi, E. R. Enzymatic conversion of D-glucose to D-fructose. *Science* **125**, 648–649 (1957).
- Bhosale, S. H., Rao, M. B. & Deshpande, V. V. Molecular and industrial aspects of glucose isomerase. *Microbiol. Rev.* **60**, 280–300 (1996).
- Adams, M. J. et al. Structure of rhombohedral 2 zinc insulin crystals. *Nature* **224**, 491–495 (1969).
- Vijayan, M. The story of insulin crystallography. *Curr. Sci.* **83**, 1598–1607 (2002).

45. Banting, F. G. & Best, C. H. The internal secretion of the pancreas. *Indian J. Med. Res.* **125**, L251 (2007).
46. Spiliopoulou, M. et al. Applications of X-ray powder diffraction in protein crystallography and drug screening. *Crystals* **10**, 54 (2020).
47. Brange, J. & Langkjær, L. *Insulin Structure and Stability*, 315–350 (Springer US, 1993).
48. Triandafillidis, D. P. et al. The T₂ structure of polycrystalline cubic human insulin. *Acta Crystallogr. Sect. D. Struct. Biol.* **79**, 374–386 (2023).
49. Gursky, O., Li, Y., Badger, J. & Caspar, D. Monovalent cation binding to cubic insulin crystals. *Biophys. J.* **61**, 604–611 (1992).
50. Gursky, O., Fantano, E., Bhyravbhata, B. & Caspar, D. Stereospecific dihaloalkane binding in a pH-sensitive cavity in cubic insulin crystals. *Proc. Natl Acad. Sci.* **91**, 12388–12392 (1994).
51. Murray, E., Draper, L. A., Ross, R. P. & Hill, C. The advantages and challenges of using endolysins in a clinical setting. *Viruses* **13**, 680–701 (2021).
52. Liu, L., Baase, W. A. & Matthews, B. W. Halogenated benzenes bound within a non-polar cavity in T4 lysozyme provide examples of I-S and I-Se halogen-bonding. *J. Mol. Biol.* **385**, 595–605 (2009).
53. Eriksson, A., Baase, W., Wozniak, J. & Matthews, B. A cavity-containing mutant of T4 lysozyme is stabilized by buried benzene. *Nature* **355**, 371–373 (1992).
54. Bouvignies, G. et al. Solution structure of a minor and transiently formed state of a T4 lysozyme mutant. *Nature* **477**, 111–114 (2011).
55. Mondal, J., Ahalawat, N., Pandit, S., Kay, L. E. & Vallurupalli, P. Atomic resolution mechanism of ligand binding to a solvent-inaccessible cavity in T4 lysozyme. *PLoS Comput. Biol.* **14**, e1006180 (2018).
56. Bar-Even, A. et al. The moderately efficient enzyme: evolutionary and physicochemical trends shaping enzyme parameters. *Biochemistry* **50**, 4402–4410 (2011).
57. Carslaw, H. & Jaeger, J. *Conduction of Heat in Solids* (Oxford University Press, 1959).
58. Krinski, S., Kazmierczak, M. & Thorne, R. Heat transfer from protein crystals: implications for flash-cooling and x-ray beam heating. *Biol. Crystallogr.* **59**, 697–708 (2003).
59. Martin, R. W. & Zilm, K. W. Preparation of protein nanocrystals and their characterization by solid state NMR. *J. Magn. Reson.* **165**, 162–174 (2003).
60. Diao, J. Crystallographic titration of cubic insulin crystals: pH affects 13 switching and sulfate binding. *Acta Crystallogr. Sect. D Biol. Crystallogr.* **59**, 670–676 (2003).
61. Vallurupalli, P., Hansen, D. F., Lundström, P. & Kay, L. E. CPMG relaxation dispersion NMR experiments measuring glycine 1H and 13C chemical shifts in the 'invisible' excited states of proteins. *J. Biomol. NMR* **45**, 45–55 (2009).
62. Vonrhein, C. et al. Data processing and analysis with the autoPROC toolbox. *Acta Crystallogr. Sect. D Biol. Crystallogr.* **67**, 293–302 (2011).
63. Tickle, I. J. et al. STARANISO. (Cambridge, United Kingdom: Global Phasing Ltd, 2016). <http://staraniso.globalphasing.org/cgi-bin/staraniso.cgi>.
64. Kabsch, W. xds. *Acta Crystallogr. Sect. D Biol. Crystallogr.* **66**, 125–132 (2010).
65. Agirre, J. et al. The CCP4 suite: integrative software for macromolecular crystallography. *Acta Crystallogr. Sect. D Struct. Biol.* **79**, 449–461 (2023).
66. McCoy, A. J. et al. Phaser crystallographic software. *J. Appl. Crystallogr.* **40**, 658–674 (2007).
67. Vagin, A. A. et al. REFMAC5 dictionary: organization of prior chemical knowledge and guidelines for its use. *Biol. Crystallogr.* **60**, 2184–2195 (2004).
68. Adams, P. D. et al. Recent developments in the PHENIX software for automated crystallographic structure determination. *J. Synchrotron Radiat.* **11**, 53–55 (2004).
69. Emsley, P., Lohkamp, B., Scott, W. G. & Cowtan, K. Features and development of Coot. *Acta Crystallogr. Sect. D Biol. Crystallogr.* **66**, 486–501 (2010).
70. Schrodinger, L. The PyMOL molecular graphics system. *Version 1*, 8 (2015).
71. Liebschner, D. et al. Polder maps: improving OMIT maps by excluding bulk solvent. *Acta Crystallogr. Sect. D Struct. Biol.* **73**, 148–157 (2017).
72. Lang, P. T., Holton, J. M., Fraser, J. S. & Alber, T. Protein structural ensembles are revealed by redefining X-ray electron density noise. *Proc. Natl Acad. Sci. USA* **111**, 237–242 (2014).
73. Desai, P. Thermodynamic properties of nickel. *Int. J. Thermophys.* **8**, 763–780 (1987).
74. Wagner, W. & Prüss, A. The IAPWS formulation 1995 for the thermodynamic properties of ordinary water substance for general and scientific use. *J. Phys. Chem. Ref. Data* **31**, 387–535 (2002).
75. Mitchell, H. M. et al. Process modelling of protein crystallisation: a case study of lysozyme. *Chem. Eng. Res. Des.* **192**, 268–279 (2023).
76. Kollie, T. Measurement of the thermal-expansion coefficient of nickel from 300 to 1000 K and determination of the power-law constants near the Curie temperature. *Phys. Rev. B* **16**, 4872 (1977).
77. Leung, A. K., Park, M. M. & Borhani, D. W. An improved method for protein crystal density measurements. *Appl. Crystallogr.* **32**, 1006–1009 (1999).
78. White, E., Tan, W. H., Ang, J. M., Tait, S. & Litster, J. The density of a protein crystal. *Powder Technol.* **179**, 55–58 (2007).
79. Powell, R., Tye, R. & Hickman, M. The thermal conductivity of nickel. *Int. J. Heat. Mass Transf.* **8**, 679–688 (1965).
80. Fujiwara, S., Maki, S., Maekawa, R., Tanaka, S. & Hagiwara, M. Measurements of thermal conductivity and thermal diffusivity of hen egg-white lysozyme crystals and its solution using the transient short hot wire method. *Int. J. Thermophys.* **38**, 123 (2017).

Acknowledgements

The authors gratefully acknowledge the support provided by the Max Planck Society. We would like to thank our colleague A.R. Pearson for critical reading of the manuscript. PM acknowledges support from the Deutsche Forschungsgemeinschaft (DFG) via grant No. 451079909 and from a Joachim Herz Stiftung add-on fellowship. ES acknowledges support by the DFG via grant No. 458246365, and by the Federal Ministry of Education and Research, Germany, under grant number 01KI2114. Funded by the European Union (ERC, DynaPLIX, SyG-2022 101071843). Views and opinions expressed are, however, those of the author(s) only and do not necessarily reflect those of the European Union or the European Research Council Executive Agency (ERCEA). Neither the European Union nor the granting authority can be held responsible for them.

Author contributions

M.S. performed protein expression, purification, prepared the protein crystals, performed data collection, analyzed the data and wrote the manuscript; C.E.H. performed protein expression, purification, crystallization, data collection and analysis; M.K., J.-P.L., H.S. and F.T. designed and validated the *spitrobot-2* and the humidity flow device (HFD) and designed the LabView interfaces; F.T., P.M. and E.C.S. designed the experiment and conceptualized the *spitrobot-2*; E.C.S. supervised and wrote the manuscript; all authors discussed and corrected the manuscript.

Funding

Open Access funding enabled and organized by Projekt DEAL.

Competing interests

On 10 March 2022, a patent application has been filed under the number EP22161384. The remaining authors declare no competing interests.

Additional information

Supplementary information The online version contains supplementary material available at

<https://doi.org/10.1038/s42004-025-01784-9>.

Correspondence and requests for materials should be addressed to Friedjof Tellkamp, Pedram Mehrabi or Eike C. Schulz.

Peer review information *Communications Chemistry* thanks Tom Caradoc-Davies, Rosemary Young, Antoine Royant, and Nicolas Caramello for their contribution to the peer review of this work.

Reprints and permissions information is available at <http://www.nature.com/reprints>

Publisher's note Springer Nature remains neutral with regard to jurisdictional claims in published maps and institutional affiliations.

Open Access This article is licensed under a Creative Commons Attribution 4.0 International License, which permits use, sharing, adaptation, distribution and reproduction in any medium or format, as long as you give appropriate credit to the original author(s) and the source, provide a link to the Creative Commons licence, and indicate if changes were made. The images or other third party material in this article are included in the article's Creative Commons licence, unless indicated otherwise in a credit line to the material. If material is not included in the article's Creative Commons licence and your intended use is not permitted by statutory regulation or exceeds the permitted use, you will need to obtain permission directly from the copyright holder. To view a copy of this licence, visit <http://creativecommons.org/licenses/by/4.0/>.

© The Author(s) 2025